

dental



clean and quiet

air

JUN-AIR®

clean air for your dental clinic




JUN-AIR introduced the first compressor almost 50 years ago and started the development of the technology and design, which has made JUN-AIR compressors the number one choice in dental clinics worldwide.

JUN-AIR invests in the latest technology. Consequently, the compressors have a compact and elegant design, reflecting quality and they are, therefore, highly suitable for installation at the place of use.

Performance and design will always play a key role in the further development of JUN-AIR compressors. At the same time, concern for the environment, low energy consumption, minimum maintenance and user-friendly operation are given high priority.

JUN-AIR supplies clean and quiet air - a complete compressed air solution for your dental clinic.

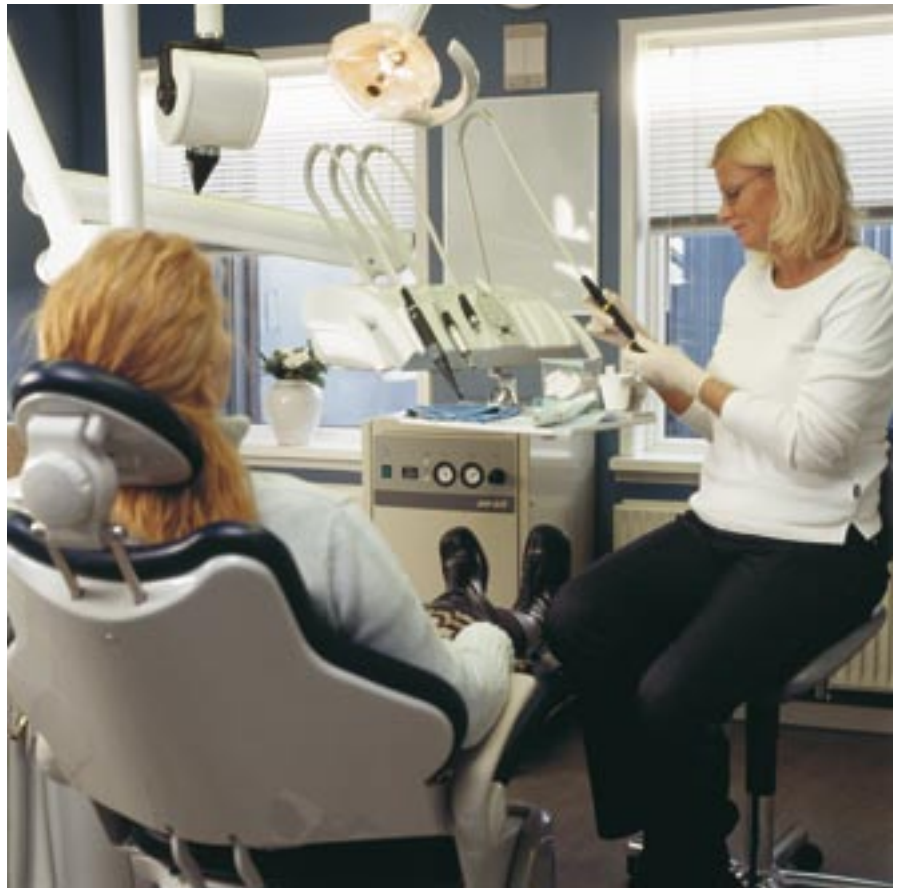

Sophisticated use of clean and dry compressed air within the dental industry results in an increasing demand for high quality air compressors.

JUN-AIR's extensive range complies with the company's traditional values - low noise level, reliability and long lifetime.

The standard products range from small models suitable for single users to large models suitable for clinics with up to 10 users.

In addition to the standard range, JUN-AIR offers a variety of accessories and large customized units, meeting the requirements of the individual customer and making JUN-AIR compressors the ideal choice for every dental clinic.

a powerful range



Model OF1202-40BD3, ideal for clinics with up to 5 users



Model OF302-25B



Model 2xOF302-40MD2

Oil-less compressors

All manufacturers of air-powered dental equipment (e.g. dental units, handpieces, micromotors, autoclaves, other hygiene products, equipment for dental laboratories, etc.), prescribe the use of clean air free from oil, particles and humidity.


The oil-less JUN-AIR compressors meet these requirements. The compressors are quiet in operation and are rated for continuous duty. They are virtually maintenance-free with a long life expectancy. The compressor receivers are internally coated to resist corrosion and the compressors are suitable for single and multiple surgeries. Moreover, the units are available with adsorption air dryer to obtain the ultimate air quality.

Oil-less in cabinet

In cases where the compressor needs to be located inside the clinic, absolute quietness is a necessity. The optimum dental solution in such cases is the M series where a metal cabinet reduces the noise by 75%, resulting in a noise level of only 47-60 dB(A).


The cabinet versions are equipped with rubber castors for easy relocation and service.

Please consult your local JUN-AIR distributor for more information on the M series.

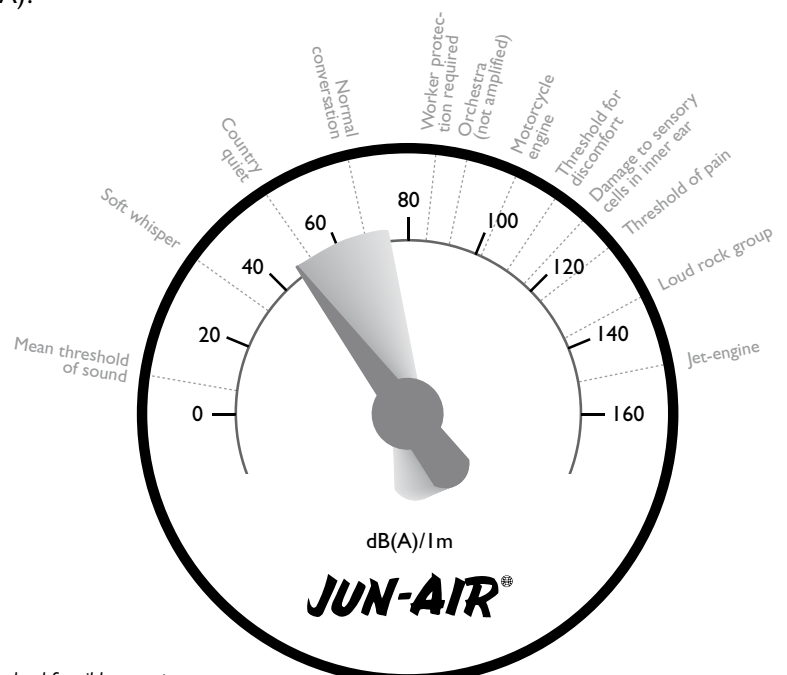

Noise level for oil-less compressors


clean and dry air

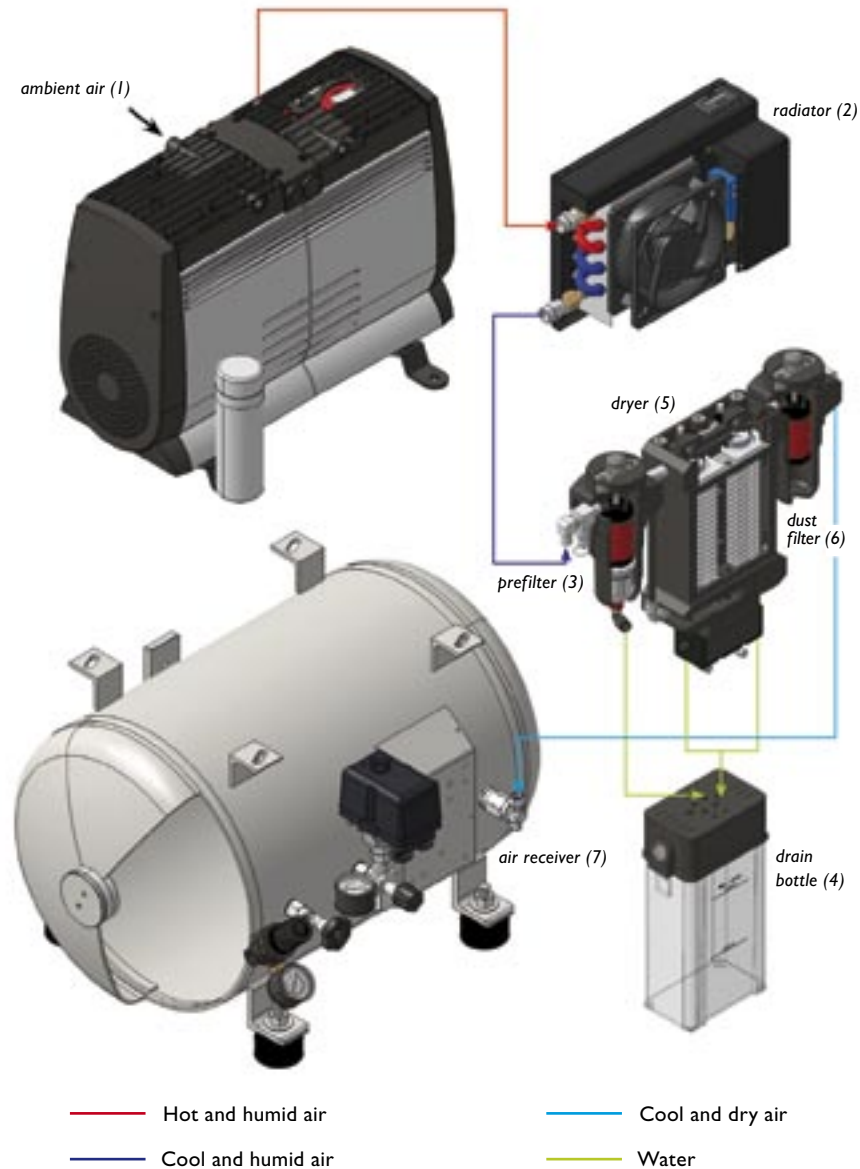

After compression (1), the air passes a radiator (2) to force as much vapour as possible to condense. The condensate (water) is trapped in the 0.01µm prefilter (3) and led to a bottle (4) for collection. The air enters the double-column adsorption air dryer (5), where the remaining vapour is removed. After having passed the dryer, the air is stored in the receiver (7), hygienically coated inside.



Sample from a compressor without dryer - moisture in the air provides ideal conditions for bacterial growth.

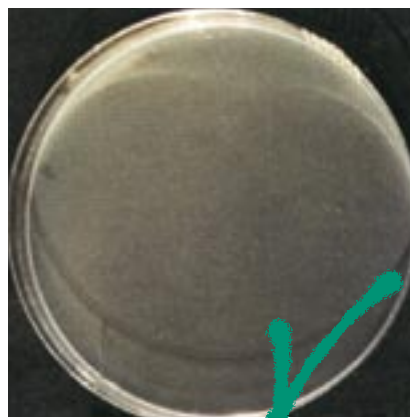

Sample taken from an oil-less compressor with dryer - as the air is dry, the risk of bacterial growth is eliminated.

Atmospheric air contains large amounts of water vapour. Unless the compressor is supplied with an efficient dryer, the major part of this water vapour will condense and change from vapour into liquid drops. This water together with the heat generated during compression, constitutes ideal conditions for the growth of microorganisms in the compressed air. Such bacteria are likely to be transferred to patients through handpieces and/or syringes.

In addition to the risk of causing contamination, water, dirt and oil in the compressed air also harm the adhesion of filling materials. And if the water is removed from the compressed air, the risk of corrosion of the equipment is eliminated.

The double-column adsorption air dryer immediately removes the water vapour, and conditions for survival of bacteria and development of corrosion are eliminated.



The dry air has a dewpoint of -40°C, which means that the air must be cooled to this temperature to condensate.

The compressed air stored for use is always 100% dry, clean and free from bacteria.

The compressor and the adsorption air dryer are designed for continuous duty and do not require standstill for regeneration. No excess capacity is required when designing an efficient compressed air system for the dental clinic based on JUN-AIR technology. Absolutely clean and dry air is always available.

When equipped with an adsorption air dryer, the JUN-AIR compressor provides compressed air, which is in accordance with the requirements of Pharmacopeia; the standard for medical air.

Basic compressors suitable for 1 dental unit

Model	OF302-25B	OF302-25BD2		
				
Number of dental units	1	1		
Mounted with dryer	No	Yes		
Voltage	V	230	230	
Frequency	Hz	50	50	
Flow:				
(FAD) @ 0 bar	l/min	108	-	
(FAD) @ 5 bar	l/min	48	38 ²⁾	
(FAD) @ 8 bar ¹⁾	l/min	38	30 ²⁾	
Tank size	litres	25	25	
Weight	kg	28	33	
Noise level	dB(A)/1m	65	65	
Dimensions (l x w x h)	mm	380 x 380 x 610	510 x 460 x 610	

Basic compressors suitable for 2-3 dental units

Model	OF1201-25HB	OF1201-25HBD2	2xOF302-40B	2xOF302-40BD2
				
Number of dental units	2	2	3	3
Mounted with dryer	No	Yes	No	Yes
Voltage	V	230	230	230
Frequency	Hz	50	50	50
Flow:				
(FAD) @ 0 bar	l/min	146	-	216
(FAD) @ 5 bar	l/min	87	69 ³⁾	95
(FAD) @ 8 bar ¹⁾	l/min	65	52 ³⁾	76
Tank size	litres	25	40	40
Weight	kg	43	53	46
Noise level	dB(A)/1m	77	77	68
Dimensions (l x w x h)	mm	634 x 413 x 626	634 x 566 x 626	560 x 450 x 610
			560 x 450 x 610	560 x 540 x 630

¹⁾ Higher pressure available upon request

²⁾ Neutral required






³⁾ Estimated value. Min. pressure required to operate dryer: 6 bar

⁴⁾ 3-phase units are approx. 100 mm wider than 1-phase units

Basic compressors suitable for 5-10 dental units

Model	OFI202-40B	OFI202-40BD3	2xOFI202-150B	2xOFI202-150BD6
				
Number of dental units	5	5	10	10
Mounted with dryer	No	Yes	No	Yes
Voltage	V	230/3x400 ²⁾	230/3x400 ²⁾	230/3x400 ²⁾
Frequency	Hz	50	50	50
Flow:				
(FAD) @ 0 bar	l/min	290	-	-
(FAD) @ 5 bar	l/min	160	128 ³⁾	256 ³⁾
(FAD) @ 8 bar ¹⁾	l/min	130	104 ³⁾	208 ³⁾
Tank size	litres	40	40	150
Weight	kg	59	70	129
Noise level	dB(A)/1m	76	76	79
Dimensions (l x w x h)	mm	556 x 443 x 678 ⁵⁾	556 x 581 x 678 ⁵⁾	1272 x 530 x 871 ⁵⁾

Cabinet compressors

Model	OF302-25M/MD2	OFI201-25M/MD2	2xOF302-40M/MD2	OFI202-40M/MD3	2xOFI202-40M/MD6
					
Number of dental units	1	2	3	5	8-10
Mounted with dryer	No Yes	No Yes	No Yes	No Yes	No Yes
Voltage	V	230 230	230 230	230/3x400 ²⁾ 230/3x400 ²⁾	230/3x400 ²⁾ 230/3x400 ²⁾
Frequency	Hz	50 50	50 50	50 50	50 50
Flow:					
(FAD) @ 0 bar	l/min	108 -	146	216 -	290 -
(FAD) @ 5 bar	l/min	48 38 ³⁾	87 69 ³⁾	96 77 ³⁾	160 128 ³⁾
(FAD) @ 8 bar ¹⁾	l/min	38 30 ³⁾	65 52 ³⁾	76 61 ³⁾	130 104 ³⁾
Tank size	litres	25 25	25 25	40 40	40 40
Weight	kg	78 83	88 98	107 112	116 121
Noise level	dB(A)/1m	47 47	58 58	50 50	60 60
Dimensions (l x w x h)	mm	720 x 460 x 860	720 x 460 x 860	670 x 650 x 860	670 x 650 x 860

¹⁾ Higher pressure available upon request

²⁾ Neutral required

³⁾ Estimated value. Min. pressure required to operate dryer: 6 bar

⁵⁾ 3-phase units are approx. 100 mm wider than 1-phase units

Denmark

JUN-AIR International A/S
Sundsholmen 3-5
Postboks 214
9400 Nørresundby
Denmark

Phone: +45 96 32 36 00
Telefax: +45 96 32 36 01
E-mail: info@jun-air.dk
Internet: www.jun-air.com

USA

JUN-AIR USA Inc.
1020 Noel Avenue, Suite A
Wheeling
IL 60090
USA

Phone: +1 847 215 9444
Telefax: +1 847 215 9449
E-mail: info@jun-air.com
Internet: www.jun-air.com

The Netherlands

JUN-AIR Benelux B.V.
Gelderlandhaven 5 B
3433 PG Nieuwegein
The Netherlands

Phone: +31 (0)30 608 3010
Telefax: +31 (0)30 608 3015
E-mail: info@jun-air.nl
Internet: www.jun-air.nl

United Kingdom

GAST GROUP Ltd
Unit 11, The I O Centre
Nash Road
Redditch, B98 7AS
United Kingdom

Phone: +44 (0)1527 50 4040
Telefax: +44 (0)1527 52 5262
E-mail: gastgroup.uk@idexcorp.com
Internet: www.gastmfg.com
www.jun-air.com

France

JUN-AIR France S.A.S.
Village Entreprises Saône Mont d'Or
444, rue des Jonchères
69730 Genay
France

Phone: +33 (0)4 37 40 82 70
Telefax: +33 (0)4 37 40 82 79
E-mail: info@jun-air.fr
Internet: www.jun-air.fr

Germany

JUN-AIR Deutschland GmbH
Kornkamp 16
22926 Ahrensburg
Germany

Phone: +49 4102 4953 0
Telefax: +49 4102 4953 45
E-mail: info@jun-air.de
Internet: www.jun-air.de

