



Softly

Mixer for pre-measured amalgam capsules

Softly is a universal high speed mixer for pre-measured capsules (amalgam, glass ionomer). It is controlled by a microprocessor which ensures precise repeatability of the selected time to ensure you a controlled consistent mix.

The aesthetically-pleasing ergonomic design of the **Softly** allows for easy use and selection of the timer display.

The width of the oscillating range (25 mm), and the mixing frequency (4000 vibrations/min. \pm 50 per min.) and your selection of the correct time required, provide perfect homogeneity and an ideal consistency for the type of amalgam that you choose.


Softly is very stable, extremely quiet, and very easy to clean and maintain.

To guarantee totally safe operation, **Softly** is deactivated when you open the front cover to access the capsule.


Softly is entirely made of injection-moulded ABS material. It is in compliance with all electrical and mechanical safety standards.



Precise repeatability of the selected time

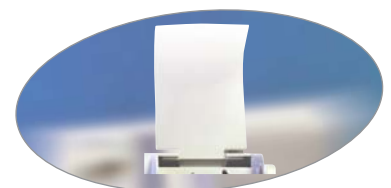

Disactivated when the front cover is open

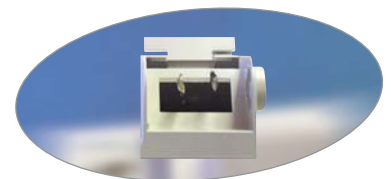

Ideal consistency for the type of filling material you choose

Advantages

- High mixing frequency
- Noiseless operation
- Excellent stability
- Full electrical safety
- Ergonomic design
- Maintenance-free
- Easy to clean



Technical specifications

Nominal voltage	.230 V- /50 Hz - 110 V- /60Hz upon request
Absorbed power	.50W
Fuse	.500 mA (230V-) - 1 AT (110 V-)
Mixing frequency	.4000 vibrations/min.
Amplitude	.25 mm
Duty	.intermittent
Programmable mixing time	.0 to 30 sec.
Dimensions	.210 x 205 x 120mm
Weight	.3.1 kg

Manufactured in compliance with the currently applicable regulations (medical devices and standards (Medical Devices EC Guideline 93/42/EEC).

